

K ONE

Design for children



Scientific operating sequence
Excellent cutting efficiency

K ONE



DESCRIPTION


Kone is a brand new type of tooth file, consisting of a blade and a handle. The front end of the blade is spiral-shaped, and there is a connecting section between the blade and the handle.

This set of files features a brand new structural design, measuring 15mm in length, suitable for root canal treatment of children's teeth. Additionally, it is compatible with the oral cavity and root canal length of children.


It can help prevent accidental damage to the child's periapical tissues during surgery and also help prevent the accidental shortening of the root canal file, which could lead to the child accidentally ingesting or inhaling it.

FEATURES

Material	NiTi
Length	15 mm
Sterilization	≤ 134°C
Package	4 pcs per box/blister assorted/single size

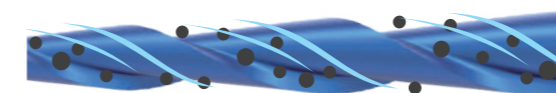

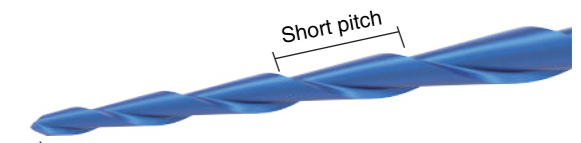

Long tail pitch design, easy to discharge chips



No cutting tip

size	Length	Torque	Speed	Color
#25/.06	15(mm)	1.50-3.0N/cm	250~350rpm	●

 Special tip design to avoid irritation and damage to the periapical tissue of primary teeth.

 The length is only 15 mm, which is suitable for most children's deciduous teeth root canal length.