

ZR

 perfect®

wireless endo motor



TOUCH

ZR TOUCH

industrially modern

exquisitely made

wireless endo motor

FEATURE

Wireless endo motor

“ZR Touch” is specially designed for dentists and used to drive root canal instruments to rotate and reciprocate. It is clinically used in root canal shaping and cleaning in root canal treatment.

	MODE	FUNCTION
M1–M5	Continuous rotation	Rotate in one direction (clockwise / counterclockwise) Speed range: 150–1000 rpm, step value 50rpm Torque range: 0.6–5.0 N·cm, step value 0.1N·cm
M6–M8	Reciprocating rotation	Speed range: 300rpm Angle range: 150° forward, 30° reversal 170° forward, 50° reversal 50° forward, 170° reversal



Adjustable digital LCD display

The ZR TOUCH has an adjustable LCD screen for both left and right hand use.




Preset file program

Built-in preset programs for different brands and models make it more convenient to use.



360° adjustable contra-angle

360° adjustable reverse-angle design allows your operation more flexible and simple.

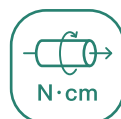

Large torque

Featured with larger torque.
The maximal torque can come to 5N·cm.