


T FLEX

 *perfect*

Innovative / Efficient

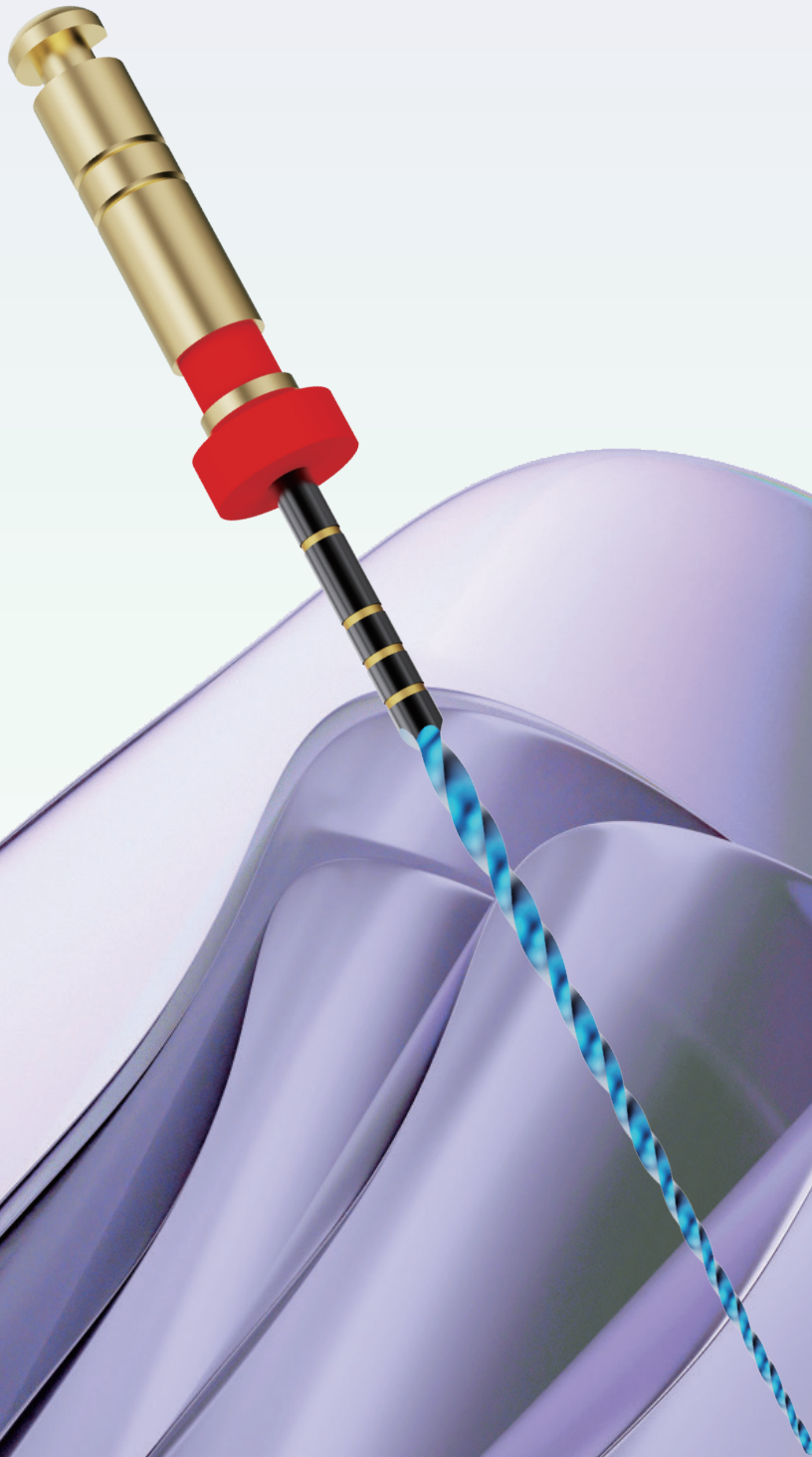

T FLEX



Innovative / Efficient



Size	Length	Torque
20/ .10	17mm	2.5N·CM
15/ .03	21,25,31mm	1.5N·CM
20/ .04	21,25,31mm	2.0N·CM
25/ .04	21,25,31mm	2.5N·CM
25/ .06	21,25,31mm	2.5N·CM

Cross section:

Rotary speed: 300-400RPM

Advantages

Extraordinary design in Flexibility

T-Flex has much more cutting efficiency and flexibility compared with ordinary rotary files.

T-Flex files are made of high purity black nickel titanium alloy which has low content of carbon and oxygen impurities.

After files heat treatment process, T-Flex conducts rotary cutting surrounding the curved axis of canals which is better adapt to the inherent anatomical structure of canals.

Perfectly for 99% root canals

T-Flex can be perfectly suitable for simple and complicated of root canals, better for curvature.

T-Flex, with small resilience can greatly avoid forming ledges, central deviation, transportation and deformation by cutting the inner and outer walls of root canals evenly.