

GP

*Precise temperature control •
Smooth extrusion of gutta-percha •*




FILL


- Supports 360° rotation
- Long - lasting battery life

• Smooth gutta - percha extrusion without clogging

The precise internal structure and specially designed internal channels enable the gutta - percha to be smoothly extruded without jamming.

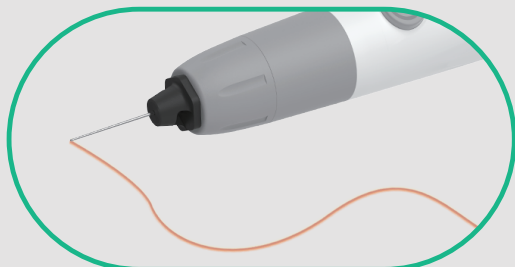

Pre-bendable needle tip

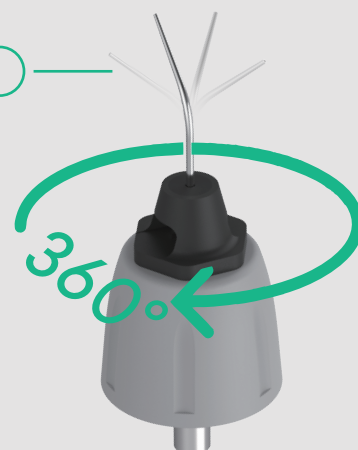

Use it in combination with Obturator Pellets



• High - temperature sterilization

The glue - injection needle is detachable
It can be sterilized by high - temperature and high - pressure steam Reusable

• Two adjustable speeds



Normal speed



Accelerated gutta-percha extrusion

• Four adjustable temperature settings

140°C


160°C

180°C

200°C

• Real-time display of the gutta-percha quantity

Equipped with an OLED display screen
Graphically display the remaining quantity of gutta-percha
Remaining quantity of gutta-percha

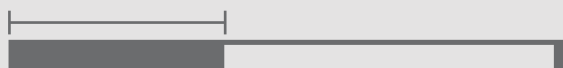

• Large - capacity battery with long - lasting endurance

Equipped with a 2600mAh high - capacity battery
Reduces the number of charging times and improves work efficiency

Model	GP Fill
Battery capacity of the handle	3.7V , 2600mAh
Power adapter	Type-c , 5V 2A
Temperature range	140°C, 160°C, 180°C, 200°C
Charging time	3h

