






RT

SCRAPER



-  Efficiently remove
-  Replaceable working tips
-  High-quality materials
-  High safety

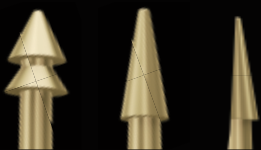

Advantages

Efficiently remove

With a variety of working tip designs, precisely penetrate deep into the root canal to remove the apical gutta-percha, clean up the residual pulp, debris, etc.

Replaceable working tips

Select and replace the appropriate working tips according to the diameter of the root canal and the filling condition of the gutta-percha.

High-quality materials


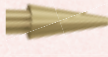
Excellent corrosion resistance, antibacterial properties and strength. This ensures stable performance during multiple uses and disinfection processes, and it is not prone to deformation.


High safety

No need to frequently adjust the angle when removing foreign objects. It will not damage the root canal wall.



Informations

Tip shape	Model	Tip diameter	Length
 Round lance shape (360°)	0.5S	0.50mm	21mm
 Round lance shape (360°)	0.3S	0.34mm	21mm
 Semicircular lance shape (180°)	0.3S-CL	0.34mm	21mm
	0.3S-CR	0.34mm	21mm
	0.3L-CL	0.34mm	28mm
	0.3L-CR	0.34mm	28mm

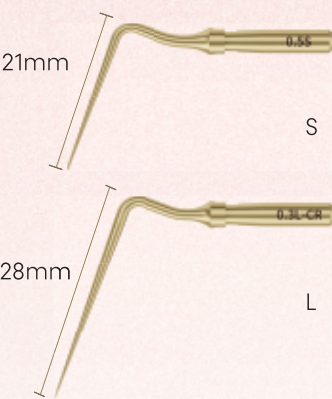

Instructions

0.3S	0.5S	0.3S-CR/0.3S-CL	0.3L-CR/0.3L-CL
			
Narrower root canals	Wider root canals.	Small diameters and areas with irregular shapes close to the canal walls.	Root canals that are long in length, small in diameter, and areas with irregular shapes close to the canal walls.
Remove the gutta-percha in the root canal, clean up the residual pulp, debris, etc.		Remove the gutta-percha and root filling materials in the root canal, as well as the necrotic pulp tissues, dentin debris, etc. in the areas that are difficult to reach.	