



Contact us


Ms. Oscar
E-mail: sales@dental-perfect.com
WhatsApp: 0086 158 1407 8725

design for children

"Kiddo" is a brand new type of file of deciduous teeth comprises a blade and a handle. The front end of the blade is spiral-shaped, and its cross-section shape is triangular. A connecting section is arranged between the blade and the handle.

This set of file system adopts new structure design, 16 mm long, contains four 4 taper and two 6 taper file, is able to apply to children's teeth root canal treatment, at the same time, appropriate to the child's mouth and root canal length.

It can help avoiding operation in the process of accident harm to children's periapical tissue, also can help avoiding the root canal file accidental shorting ingestion aspiration in children.

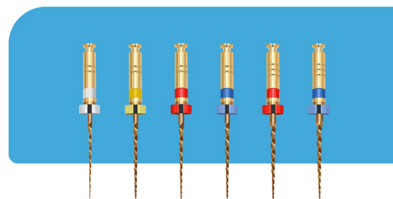

KIDDO



Special tip design to avoid irritation and damage to the periapical tissue of primary teeth

16
mm

The length is only 16mm, which is suitable for most children's deciduous teeth root canal length

Shenzhen PERFECT Medical Instruments Co., Ltd.

1F, Bldg C Plant No. 3, Jixia Zaohekeng Industrial Zone, Nanwan St., Longgang Dist., Shenzhen, China

www.dental-perfect.com

KIDDO



www.dental-perfect.com

Scientific operating sequence

Excellent cutting efficiency



Aims to be the most positive influential dental healthcare brand. To create outstanding global dental brand, promoting innovations through state-of-the-art manufacturing.

cross section



Material: Gold NiTi Package: 6 pcs per box/blister Assorted/single size

Length: 16mm Sterilization: ≤134°C

Size	Length	Torque	Taper	Speed	ISO Color
#15	16(mm)	1.50N/cm	.04	200-350(rpm)	White
#20	16(mm)	1.50N/cm	.04	200-350(rpm)	Yellow
#25	16(mm)	1.50N/cm	.04	200-350(rpm)	Red
#30	16(mm)	1.50N/cm	.04	200-350(rpm)	Blue
#25	16(mm)	1.50N/cm	.06	200-350(rpm)	Red
#30	16(mm)	1.50N/cm	.06	200-350(rpm)	Blue

