


WeApex



- Measurement & Pulp Vitality Test
- 4.3" IPS Touch Screen
- Bluetooth connection
- Wireless charging

Multifunctional integration

WeApex enables accurate root canal length measurement and pulp vitality testing, and can be seamlessly connected to the Perfect Motor to achieve simultaneous canal shaping and real-time measurement, greatly improving efficiency and accuracy.

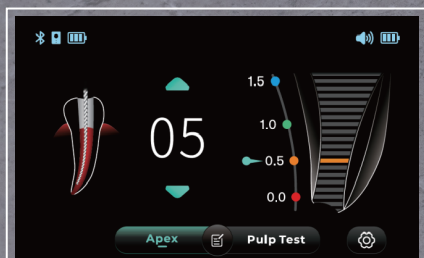

Measure the working length of the root canal.



Test pulp vitality

4.3-inch IPS screen

The 4.3-inch IPS screen provides a clear display at a glance. It features touch control, which is quick and convenient. The interface is simple and clear, enabling you to operate all functions just by touching.



Wireless charging

Apex Locator supports wireless charging. Simply place it on the Apex Dock, and it will start charging automatically.

Large battery capacity with strong endurance
Battery capacity of the base: 2600mAh. It can charge Apex Locator at least 10 times.

Precise measurement

Upgraded to the fifth-generation multi-frequency root canal measurement technology, ensuring more stable and accurate measurement.

Model	WeApex
Dimensions	43.5*21.2*14.1mm
Battery Capacity of Apex Dock	3.7V, 2600mAh
Battery Capacity of Apex Locator	3.7V, 160mAh
Power Adapter Output	Type C, 5V 2A
Screen	4.3-inch touch-control color screen
Weight	About 14g